



PETERMANN®

Technology for people of limited mobility



*"Petermann –
because every
movement counts"*

blue
MOTION
TRANSFER AIDS

alpha
TRANSFER AIDS

ONLYU
TRANSFER AIDS

Index

Your quick overview

	General
	About us 6
	Product overview8
	YouTube Channel14
	Alpha®
	Sliding Mats 19
	Transfer Sheets 22
	Raising and Lifting Aids 29
	Oval Boards 36
	Sliding Boards..... 37
	Sliding Mats Wheelchair 38
	Aids for Daily Care 41
	X-Ray Aids 44
	Rollboards 46
	OnlyU®
	Rollboards 52
	Disposable Sliding Foil 55
	Wall Rack for Foils 57
	Blue Motion®
	Slides 60
	Evacuation systems
	Evacuation System 67



Innovations

We know what you need

Petermann product ideas emerge exactly where they will be required and used later on – in day-to-day life. Every Petermann product has been developed to meet the requirements of caregivers and patients who need to master a sequence of movements. The Petermann teams are involved in the daily activities of caregivers and see every day the heavy the daily burden and challenge they face when caring for people with limited mobility.

Quality

Safe products for proficient care –
Made in Germany

Many years of experience, the latest production technologies, innovative materials and the highest degree of care in production – These are some of the criteria that highlight the extraordinary quality of the Petermann aids.

Continual striving for perfection, finding even better materials and processing techniques – linked to systematic quality tests carried out in both our own laboratories and also in independent testing centres – enable Petermann GmbH to guarantee that all products fully match requirements for medical products. They are therefore free of any hazardous substances and also withstand the toughest stresses and demands.



Petermann products have above average life cycles. With a period of up to 20 years, the guarantees go well beyond the legal stipulations. Most products are also suitable for heavy loads, with some products able to withstand loads of up to 400 kg.

Cost effectiveness

The best offer for you

The modern health economy regards cost effectiveness as the sum of many different factors. Petermann GmbH offers you fair purchase prices for aids with long service lives and also for patient-related disposable products plus many options for taking advantage of Petermann's expertise product handling and training.

Petermann GmbH is familiar with the requirements in the care sector and it knows that products must also continue to guarantee security of care after many washes and disinfection processes for both users and patients.

Back-friendly movement for patients and relatives

Allow yourself to be "helped"
and pay attention to your back!

Petermann products helps all caregivers carry out their activities in a back-friendly manner while gently moving and transferring patients. *"It is necessary to make things as simple as possible. But no simpler than that!"* reads the Petermann motto; this means making the products available to caregivers and relatives providing care with very different levels of training and experience both simply and effectively.

In order to ensure that everyone can work with our products as efficiently as possible, we offer product training and intensive training plus complete Bella Reha one-day seminars which convey knowledge about practical handling of Petermann products and provide lots of effective, easy-to-apply advice.

Qualified expertise

Put your faith in
Petermann quality

Petermann GmbH sees itself as Germany's pioneer in developing back-friendly aids for caregivers. Experienced product specialists at Petermann GmbH have carried out many thousands of training sessions and seminars in care homes, hospitals, outpatient care facilities and medical institutions over the years. Our key mission is to bring together this knowledge and experience we have acquired at Petermann GmbH and passing it on to you as an integrated patient transfer concept.

The Petermann quality management system is tested and certified by independent bodies. This is how Petermann GmbH ensures complete compliance with the highest quality and safety standards, from production to application.


PETERMANN®
Technology for people of limited mobility




Alpha® Caremaster
Open Sliding Mat



Caremaster XL	
Item no.	PM-2097
Dimensions	185 x 70 cm
Material	PES/Carbon
Caremaster M	
Item no.	PM-2095
Dimensions	175 x 60 cm
Material	PES/Carbon
Caremaster XXL	
Item no.	PM-2094
Dimensions	220 x 80 cm
Material	PES/Carbon

S. 19


Alpha® Sliding Mat

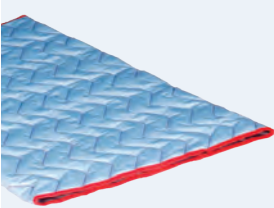

Sliding Mat	
Item no.	PM-1210
Dimensions	100 x 70 cm
Material	PES/Carbon
Sliding Mat large	
Item no.	PM-1215
Dimensions	110 x 75 cm
Material	PES/Carbon

S. 20

Alpha® Sliding Mat

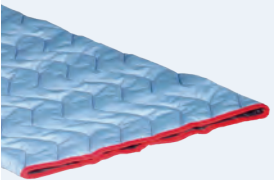

Sliding Mat	
Item no.	PM-1010
Dimensions	190 x 60 cm
Material	PES/Carbon

S. 20

Alpha® Sliding Mat



Sliding Mat	
Item no.	PM-2010
Dimensions	50 x 60 cm
Material	PES/Carbon
Sliding Mat 2-parts	
Item no.	PM-2015
Dimensions	25 x 60 cm
Material	PES/Carbon

S. 21


Alpha® Sliding Mat



Sliding Mat	
Item no.	PM-2020
Dimensions	50 x 60 cm
Material	Cotton/ natural ticking
Sliding Mat	
Item no.	PM-1220
Dimensions	100 x 70 cm
Material	Cotton/ natural ticking

S. 21

Alpha® Bed Wedges
Set 3-parts

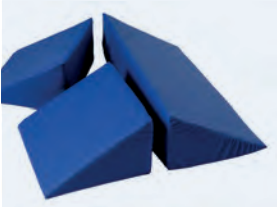

Bed Wedges	
Item no.	PM-8010
Dimensions	25 x 22 x 13 cm
2 x Bed Wedges	25 x 22 x 13 cm
Dimensions	50 x 22 x 13 cm
1 x Bed Wedges	50 x 22 x 13 cm
Material	PU-Non-Slip

S. 24

Alpha® Sliding Mat



Sliding Mat	
Item no.	PM-1420
Dimensions	140 x 90 cm
Material	Cotton/ natural ticking

S. 21

Alpha® Turn Sheet
Set 135



Turn Sheet Set 135	
Item no.	PM-1580
consisting of:	
Turn Sheet	PM-1582
Fitted Sheet	PM-1584
Non-Slip-Mat	PM-7010
Bed Wedges	PM-8010
Turn Sheet	
Item no.	PM-1582
Dimensions	150 x 200 cm
Material	Cotton/ natural ticking
Fitted Sheet	
Item no.	PM-1584
Dimensions	200 x 90 cm
Non-Slip-Mat	
Item no.	PM-7010
Dimensions	35 x 25 cm
Bed Wedges	
Item no.	PM-8010
Dimensions	25 x 22 x 13 cm
2 x Bed Wedges	25 x 22 x 13 cm
Dimensions	50 x 22 x 13 cm
1 x Bed Wedges	50 x 22 x 13 cm

S. 25

Alpha® Drehlaken
Set Air



Turn Sheet Set Air	
Item no.	PM-1699
consisting of:	
Turn Sheet	PM-2099.1
with pocket	PM-2091
Fixing Stripes-set	PM-1584
Fitted Sheet	PM-1232.1
Air Cushion	
Cover for	PM-1233
Air Cushion	PM-1234
Air Pump	PM-7010
Non-Slip-Mat	PM-8010
Bed Wedges	
Turn Sheet with pocket	
Item no.	PM-2099.1
Dimensions	200 x 100 cm
Material	Cotton/ natural ticking
Fitted Sheet	
Item no.	PM-1584
Dimensions	200 x 90 cm
Air Cushion	
Item no.	PM-1232.1
Cover for Air Cushion	
Item no.	PM-1233
Air Pump	
Item no.	PM-1234
Non-Slip-Mat	
Item no.	PM-7010
Dimensions	35 x 25 cm
Bed Wedges	
Item no.	PM-8010
2 x Bed Wedges	25 x 22 x 13 cm
1 x Bed Wedges	50 x 22 x 13 cm

S. 27

Alpha® Turn Sheet
Set Multi Turn



Turn Sheet Set Multi Turn	
Item no.	PM-1695
consisting of:	
Turn Sheet	PM-2098
Fitted Sheet	PM-1584
Non-Slip-Mat	PM-7010
Bed Wedges	PM-8010
Turn Sheet with grips	
Item no.	PM-2098
Dimensions	200 x 100 cm
Material	Cotton/ natural ticking
Fitted Sheet	
Item no.	PM-1584
Dimensions	200 x 90 cm
Non-Slip-Mat	
Item no.	PM-7010
Dimensions	35 x 25 cm
Bed Wedges	
Item no.	PM-8010
Dimensions	25 x 22 x 13 cm
2 x Bed Wedges	25 x 22 x 13 cm
Dimensions	50 x 22 x 13 cm
1 x Bed Wedges	50 x 22 x 13 cm

S. 26

Alpha® Support Belt



Support Belt XS	
Item no.	PM-6012
for waist	45–70 cm
Support Belt S	
Item no.	PM-6013
for waist	60–95 cm
Support Belt M	
Item no.	PM-6014
for waist	72–110 cm
Support Belt L	
Item no.	PM-6015
for waist	100–150 cm
Support Belt XL	
Item no.	PM-6016
for waist	140–180 cm
Material	Polyamide

S. 29

Alpha® Raising Aid /
Stand-up Support



Raising Aid / Stand-up Support with a flexible holding belt	
Item no.	PM-3015
Dimensions	60 x 20 cm
Raising Aid Standard	
Item no.	PM-3012
Dimensions	60 x 20 cm
Material	Polyamide

S. 30

Alpha® Raising Aid-Set
Magic & Leg Belt bundle



Raising Aid-Set
Item no. PM-3025
consisting of:
Raising Aid Magic PM-3015
Leg Belt PM-3005
Leg Belt
Item no. PM-3005
Dimensions 100 x 22 cm
Material PES/Carbon
Leg Belt-Extension
Item no. PM-3007
Dimensions 50 x 22 cm
Material PES/Carbon



S. 31

Alpha® Looped Belt



Looped Belt
Item no. PM-3010
Dimensions 110 x 4 cm
Material Polyester
with 2 loops including buckles

S. 30

Alpha® Lifting Mat



Lifting Mat 4 handles
Item no. PM-4012
Dimensions 110 x 45 cm
Lifting Mat 6 handles
Item no. PM-4016
Dimensions 150 x 56 cm
Material Nylon

S. 32

Alpha® Transfer Seat Swing



Transfer Seat Swing
Item no. PM-4035
Dimensions 38 x 43 x 55 cm
Material Polyamide

S. 33

Alpha® Raising Aid Quick



Raising Aid Quick
Item no. PM-6080
Dimensions 145 x 7–13 cm
Material Polyester

S. 33

Alpha® Transfer Aid Merlin



Transfer Aid Merlin, 3-parts
Item no. PMM-1000
consisting of:
Transfer Aid PMM-1100
Loop Belt PM-3010
Raising Aid PM-6160
Transfer Aid Merlin
Item no. PMM-1100
Base frame 52,9 x 81,0 cm
Footplate 41,0 x 46,5 cm
Base frame Stainless Steel

S. 35

Alpha® Support Belt Merlin



Support Belt Merlin
Item no. PM-6160
for waist 70–160 cm
Material Polyamide

S. 35

Alpha® Flexi Board Oval



Flexi Board Oval S (pair)
Item no. PM-5045
Dimensions 2 x 46 x 20 cm
Material Polyethylene
Flexi Board Oval L (pair)
Item no. PM-5040
Dimensions 2 x 85 x 40 cm
Material Polyethylene

S. 36

Alpha® Flexible Sliding Board



Sliding Board with rounded edges
Item no. PM-5030
Dimensions 74 x 32 cm
Material Polyethylene
Sliding Board
Item no. PM-5050
Dimensions 74 x 32 cm
Children's Sliding Board
Item no. PM-5020
Dimensions 58 x 32 cm
Material Polyethylene

S. 37

Alpha® Sliding Mat Wheelchair



Sliding Mat Wheelchair
Item no. PM-2070
Dimensions 50 x 45 cm
Material Cotton/
natural ticking

S. 38

Alpha® Sliding Mat Wheelchair



Sliding Mat Wheelchair
Item no. PM-2076
Dimensions 130 x 45 cm
Material Cotton/
natural ticking

S. 38

Alpha® Sliding Mat Wheelchair



Sliding Mat Wheelchair
Item no. PM-2078
Dimensions 135 x 43 cm
Material Cotton/
natural ticking

S. 39


Alpha® Turning Disc



Turning Disc
Item no. PM-5080
Diameter 40 cm
Material ABS coated
plastic fabric

S. 41

Alpha® Non-Slip Cushion

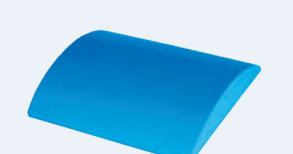

Non-Slip Cushion
Item no. PM-8050
Dimensions 40 x 50 cm
Material Polyethylene

S. 42

Alpha® Footstool



Footstool
Item no. PM-8090
Dimensions 32 x 38 x 8,5 cm
Material Polyethylene
Footstool Cover
Item no. PM-8091
Material Polyamide

S. 41

Alpha® Non-Slip-Mat



Non-Slip-Mat
Item no. PM-7010 red
Item no. PM-7010-1 blue
Dimensions 35 x 25 cm
Non-Slip-Mat large
Item no. PM-7011
Dimensions 45 x 38 cm
Material Thick film

S. 42

Alpha® Bed Ladder



Bed Ladder
Item no. PM-6090
Dimensions 170 x 8 cm
Material Polyester

S. 43


Alpha® Soft Turning Disc



Soft Turning Disc
Item no. PM-5070
Diameter 46 cm
Material Cotton/
natural ticking

S. 43


Alpha® X-Ray Mat

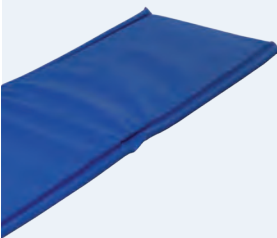

X-Ray Mat 60 cm	
Item no.	PMF-8200
Dimensions	190 x 65 cm
X-Ray Mat 85 cm	
Item no.	PMF-8220
Dimensions	190 x 85 cm
Material	PVC/PES

S. 44


Alpha® X-Ray Crawler

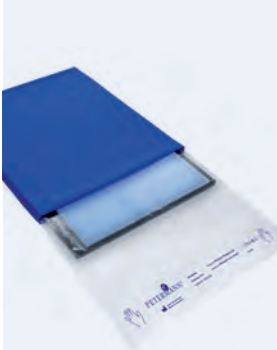

X-Ray Crawler-Set 3-parts	
Item no.	PMF-8265
2 x X-ray Crawler	PMF-8260
1 x Protection Foil	PMF-8261
X-Ray Crawler-Set 6-parts	
Item no.	PMF-8266
5 x X-ray Crawler	PMF-8260
1 x Protection Foil	PMF-8261
X-Ray Crawler	
Item no.	PMF-8260
Dimensions	41 x 53,5 cm
Material	Polyamide/PE
X-Ray Crawler-Set XL 3-parts	
Item no.	PMF-8365
2 x X-ray Crawler	PMF-8360
1 x Protection Foil	PMF-8361
X-Ray Crawler-Set XL 6-parts	
Item no.	PMF-8366
5 x X-ray Crawler	PMF-8360
1 x Protection Foil	PMF-8361
X-Ray Crawler XL	
Item no.	PMF-8360
Dimensions	47 x 53,5 cm
Material	Polyamide/PE

S. 45

Alpha® Protection Foil
Accessories for X-Ray Crawler

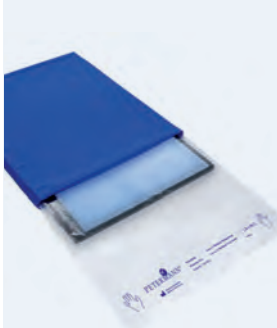

Protection Foil	
Item no.	PMF-8261
Dimensions	45 x 80 cm
Material	LDPE
Protection Foil XL	
Item no.	PMF-8361
Dimensions	51 x 80 cm
Material	LDPE
Perforated on a roll with 100 pcs	

S. 45

Alpha® Rollboard Vision



Rollboard Vision	
Item no.	PMF-9050
Dimensions	85 x 45 cm 2 x foldable
Item no.	PMF-9075
Dimensions	85 x 45 cm foldable
Item no.	PMF-9085
Dimensions	95 x 50 cm foldable
Item no.	PMF-9095
Dimensions	95 x 50 cm not foldable
Item no.	PMF-9150
Dimensions	150 x 60 cm foldable
Item no.	PMF-9155
Dimensions	150 x 60 cm not foldable
Item no.	PMF-9170
Dimensions	180 x 50 cm foldable
Item no.	PMF-9180
Dimensions	195 x 50 cm foldable
Material	Polyethylene
Rollboard	Polyamide
Cover	
Item numbers replacement covers see page 46	

S. 46


Alpha® Rollboard Vision
with cover with handles and grips



Rollboard Vision	
Item no.	PMF-9185
Dimensions	180 x 65 cm foldable with loops and handles
Replacement Cover	PMF-9186
Material	Polyethylene
Rollboard	Polyamide/PES
Cover	

S. 49

Alpha® Wall Rack
for Rollboards

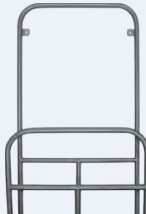

Wall Rack	
Item no.	PMF-0950
Dimensions	75 x 45 x 9 cm
Material	Stainless Steel

S. 49


OnlyU® Rollboard Vision PEM



Rollboard Vision	
Item no.	PMF-9089
Dimensions	95 x 50 cm foldable
PEM Disposable Cover	PMF-9090
Item no.	PMF-9099
Dimensions	95 x 50 cm not foldable
PEM Disposable Cover	PMF-9090
Item no.	PMF-9153
Dimensions	150 x 60 cm foldable
PEM Disposable Cover	PMF-9152
Item no.	PMF-9157
Dimensions	150 x 60 cm not foldable
PEM Disposable Cover	PMF-9152
Item no.	PMF-9175
Dimensions	180 x 50 cm foldable
PEM Disposable Cover	PMF-9172
Item no.	PMF-9189
Dimensions	195 x 50 cm foldable
PEM Disposable Cover	PMF-9172
Material Rollboard	Polyethylene
PEM Disposable Cover	
Item no.	PMF-9090
Dimensions	90 x 53 cm
Quantity	120 pcs on roll, perforated
Item no.	PMF-9152
Dimensions	150 x 61 cm
Quantity	70 pcs on roll, perforated
Item no.	PMF-9172
Dimensions	180 x 53 cm
Quantity	60 pcs on roll, perforated
Item no.	PMF-9077
Dimensions	85 x 47,5 cm
Quantity	50 pcs folded
Material	Polyethylene
Disposable Protection Foil on roll	
Item no.	PMF-0191
Dimensions	90 x 54,5 cm
Quantity	110 pcs
Item no.	PMF-0192
Dimensions	180 x 54,5 cm
Quantity	55 pcs

S. 52

OnlyU® PEM Disposable Sliding Foil

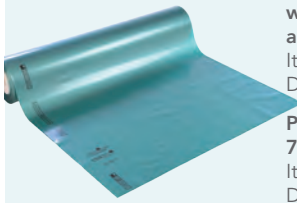

PEM Disposable Sliding Foil 60 cm	
Item no.	PM-060.60-100
Dimensions	60 cm x 100 m
PEM Disposable Sliding Foil 60 cm with dispenser box and foil cutter	
Item no.	PM-060.65-100
Dimensions	60 cm x 100 m
PEM Disposable Sliding Foil 70 cm	
Item no.	PM-060.70-100
Dimensions	70 cm x 100 m
PEM Disposable Sliding Foil 70 cm with dispenser box and foil cutter	
Item no.	PM-060.75-100
Dimensions	70 cm x 100 m
Material	Polyethylene

S. 55


OnlyU® Dispenser for Sliding Foil



Dispenser for rolls up to 70 cm	
Item no.	PM-060.30

S. 56

OnlyU® Wall Rack

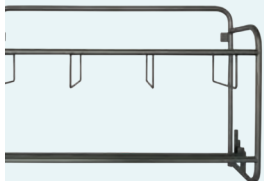

Wall Rack	
Item no.	PMF-0200.2
Dimensions	51 x 80 x 20 cm
Material	Stainless Steel
including basket for replacement foil	

S. 57

Blue Motion® Slide
in many sizes



Slide

Item no.	BM-0560*
Dimensions	52 x 60 cm
Item no.	BM-0652*
Dimensions	60 x 52 cm
Item no.	BM-0752*
Dimensions	70 x 52 cm
Item no.	BM-S-0790
Dimensions	70 x 90 cm
Item no.	BM-S-0907
Dimensions	90 x 70 cm
Item no.	BM-S-0910
Dimensions	90 x 100 cm
Item no.	BM-S-0911
Dimensions	90 x 110 cm
Item no.	BM-S-1111
Dimensions	110 x 110 cm
Item no.	BM-S-1170
Dimensions	110 x 70 cm
Item no.	BM-S-1190
Dimensions	110 x 90 cm
Item no.	BM-S-1370
Dimensions	130 x 70 cm
Item no.	BM-S-1970
Dimensions	195 x 70 cm
Material	Polyamide

*recommended for use in wet areas or in combination with the Sliding Board

S. 60

Blue Motion® Non-Slip-Mat
10 pack



Non-Slip-Mat

Item no.	BM-7005
Dimensions	25,5 x 35,5 cm
Material	Thick film

S. 62

Blue Motion® Sliding Board



Sliding Board

Item no.	BM-5000
Dimensions	74 x 24 cm
Material	Polyethylene

S. 63

Blue Motion® Sheet
with sewn-in handles



Sheet

Item no.	BM-K-1975
Dimensions	195 x 73 cm
Material	Polyamide

S. 62

Evacuation Sheet

Evacuation Sheet

Item no.	PM-EVA-7053-H
for bed size:	
Length	190–205 cm
Width	80–100 cm
Item no.	PM-EVA-7053-S
for bed size:	
Length	190–205 cm
Width	100–140 cm
Material	PVC-coated polyester fabric completely perforated

S. 66/67

Petermann YouTube Channel

Videos of our products

Our channel contains examples of how to use our repositioning and transfer aids. A range of applications for different health conditions as well as backfriendly work methods are introduced in detail.

www.pm-med.de

INFO

Please bear in mind that our videos do not replace product training!




„Take a look round
and visit our new
YouTube Channel!“



Alpha[®] Caremaster

Open Sliding Mat

The Alpha Caremaster Sliding Mat makes care work considerably easier. Shear forces can be almost completely eliminated. The Alpha Caremaster Sliding Mat is particularly versatile in its possible applications and is used for all possible position changes while lying down and for transfers in bed. It provides a secure support and an impulse-giving movement sequence without having to touch the patient in problem areas during the transfer. As the usual tunnel shape has been dispensed with here, the pull loops can be used to effortlessly move the patient to the head end and to the side without having to remove and reinsert the aid.

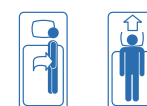

Caremaster XL

Item no.	PM-2097
Dimensions	185 x 70 cm
Material	PES/Carbon

Caremaster M

Item no.	PM-2095
Dimensions	175 x 60 cm
Material	PES/Carbon

Caremaster XXL

Item no.	PM-2094
Dimensions	220 x 80 cm
Material	PES/Carbon


Special sizes on request



Soft padding allows comfortable transfer



Head support for a gentle transfer

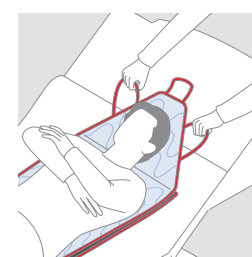


The Caremaster with three extra-strong webbed loops can also be used for very heavy persons


Alpha® Sliding Mat

Sliding Mats PES/Carbon

Alpha Sliding Sheets make it particularly easy to transfer people in need of care in bed, as they significantly reduce friction resistance. Depending on the functional profile of the person in need of care, Alpha Sliding Sheets simultaneously support the chest and pelvis for positioning and turning onto the side. Another possible application, such as repositioning a heavy person in need of care from a supine to a prone position, is possible with Alpha Sliding Mats.

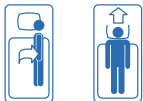

The service life of the Sliding Mat is increased by strong seams

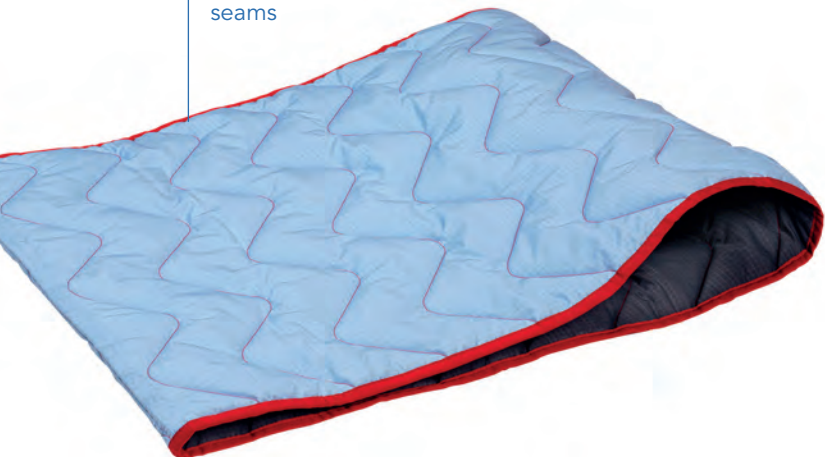

PM-2010 Sliding Mat




PM-2015 Sliding Mat 2-parts


The Alpha Sliding Mats are softly padded and also offer good protection against injuries to the heels when moving to the head end



The length of the Alpha Sliding Mat corresponds to the body size

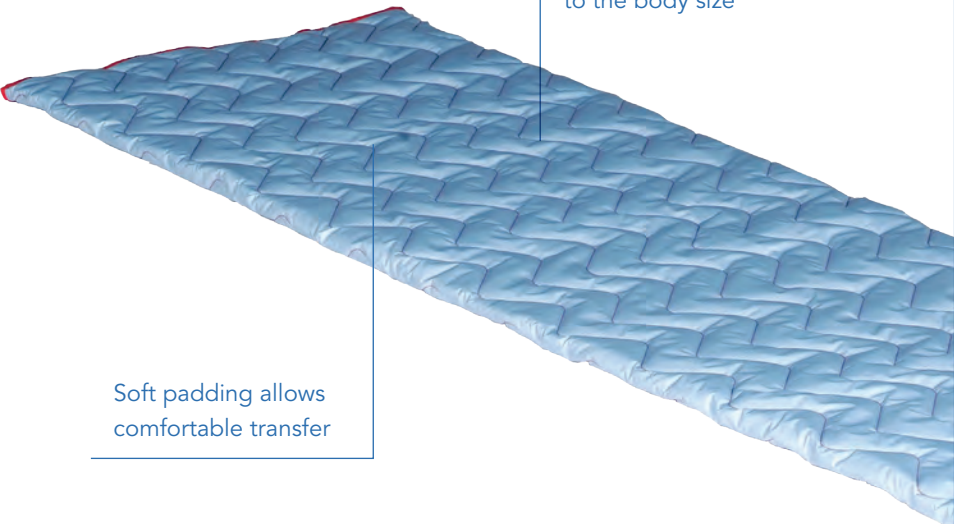



Soft padding allows comfortable transfer

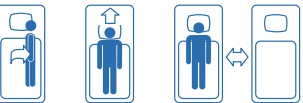

PM-1010 Sliding Mat

Alpha® Sliding Mat

Sliding Mats Cotton/natural ticking

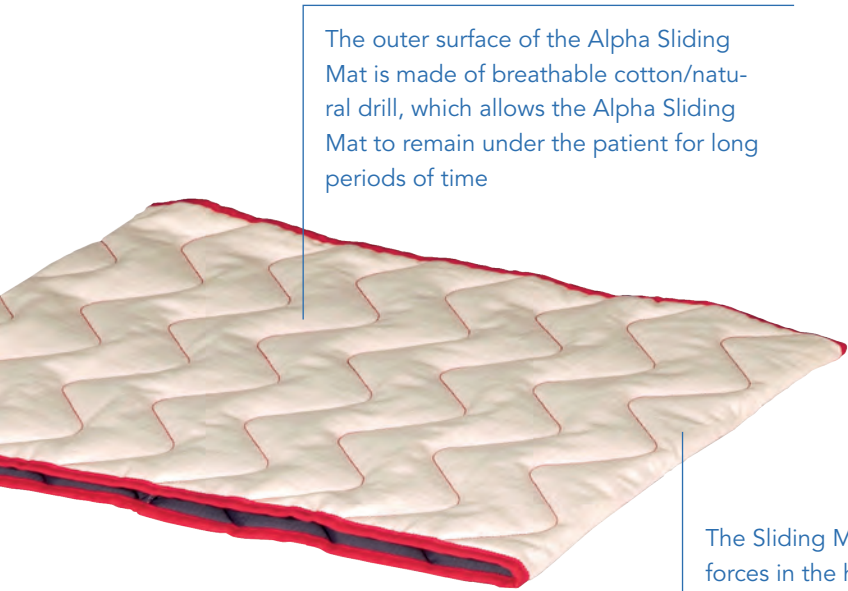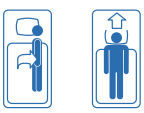

The outer surface of the Alpha Sliding Mat is made of breathable cotton/natural drill, which allows the Alpha Sliding Mat to remain under the patient for long periods of time

The Sliding Mat reduces shear and friction forces in the heel area. The Sliding Mat can be used preventively, supportively and for therapy



Item no.	PM-2020
Dimensions	50 x 60 cm
Material	Cotton/natural ticking

Sliding Mat Standard	
Item no.	PM-1220 (w/o fig)
Dimensions	100 x 70 cm
Material	Cotton/natural ticking


Alpha® Sliding Mat

Sliding Mats Cotton/natural ticking

The Alpha Sliding Mat is mainly used for turning in bed or for positioning on the side. It has four fixed straps at the top. These can be used to attach it to the slatted frame at the head of the bed if necessary. The Alpha Sliding Mat covers the entire width of the bed and 3/4 of its length.



Item no.	PM-1420
Dimensions	140 x 90 cm
Material	Cotton/natural ticking

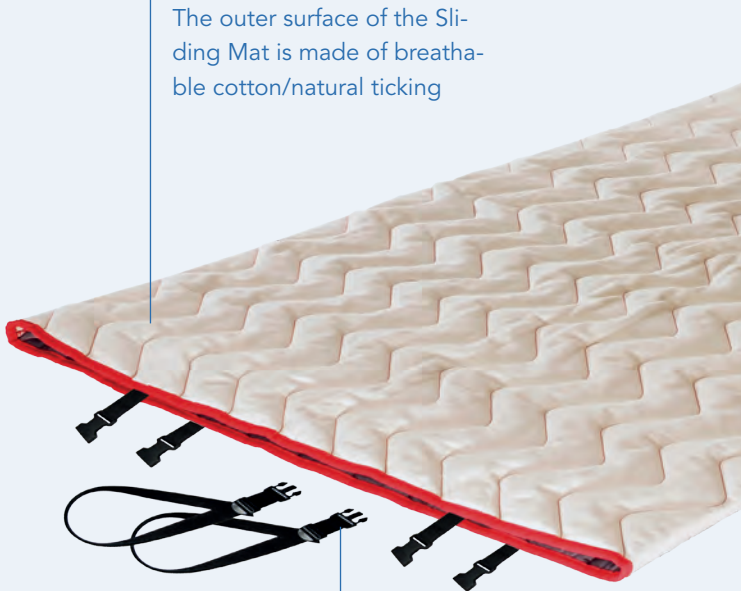

The outer surface of the Sliding Mat is made of breathable cotton/natural ticking


The Alpha Sliding Mat can be attached to the headboard of the bed with straps

Alpha® Turn Sheet Set

Turn Sheet Variants

Turn Sheet Set 135

For easy implementation of the 135° position, to relieve pressure in the back, sacrum and heel area

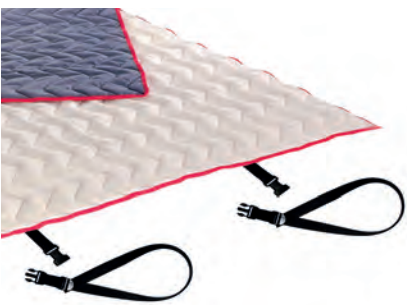

Turn Sheet Set Multi Turn

For changing position at short intervals to relieve capillary vessel pressure, as well as constant repositioning and relief of areas of the body at risk of pressure sores



Turn Sheet Set Air

Thanks to the integration of an air cushion with a variably adjustable air volume, even small movements by the nursing staff or the person in need of care themselves are sufficient to enable them to move to the side



Alpha® Turn Sheet Set

Basic equipment

The basic equipment is combined with the Turn Sheet appropriate for the respective care situation. Each basic equipment set contains:

A Fitted Sheet with a large sliding surface made of extremely low-friction material and with side stop strips that provide the necessary safety



PM-1584 Fitted Sheet

An Non-Slip-Mat that prevents your feet from sliding in bed and increases friction resistance by 100%



PM-7010 Non-Slip-Mat

Bed Wedges for active and corrective positioning of people requiring care



PM-8010 Bed Wedges

Turn Sheet Set Basic equipment

Fitted Sheet

Item no.	PM-1584
Dimensions	200 x 90 cm
Material	Cotton/natural ticking

Suitable for beds in sizes 90 x 200 cm resp. 100 x 200 cm (for a maximum mattress thickness of 22 cm)

Special sizes on request



Non-Slip-Mat

Item no.	PM-7010
Dimensions	35 x 25 cm
Material	Thick film




Alpha® Bed Wedges

Set 3-parts, part of basic equipment Turn Sheet Set

Bed Wedges were developed for the active and corrective positioning of people in need of care. Their variable use allows for targeted and purposeful positioning depending on the planned care measures.


Item no.	PM-8010
Dimensions	
2 x Bed Wedges	25 x 22 x 13 cm
1 x Bed Wedges	50 x 22 x 13 cm
Material	PU Non-Slip

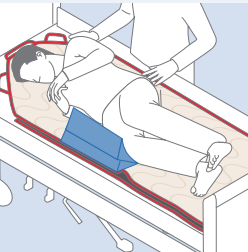




Positioning option for changing dressings on the sacrum, as well as for basic care and changing incontinence materials (secure and stable positioning up to 135°)



Positioning option for eating or to prevent slipping when sitting in the desired position in bed.

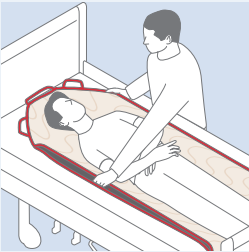

Drawings above, center, and below: Positioning options for lateral positioning and for changing position using micro-positioning, especially in combination with the Alpha Sliding Mat and Turn Sheet.

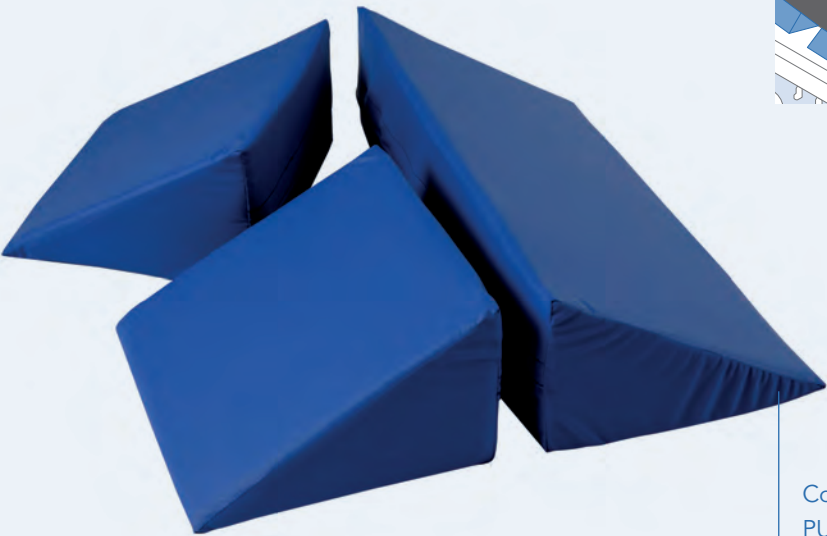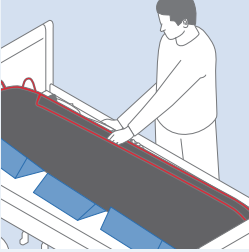
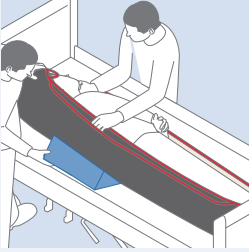


Cover made of skin-friendly PU material

Alpha® Turn Sheet



Set 135

This Turn Sheet allows a person in need of care to be moved and repositioned in bed without exerting force or causing friction. In particular, a 135° position, as an alternative to the prone position, can be easily achieved. The 135° position is particularly useful for relieving pressure on the back, sacrum, and heel areas. The use of breathable textile products means that the patient can remain on the Turn Sheet at all times. The stop strips attached to the Fitted Sheet and the Non-Slip-Mat ensure the necessary safety.

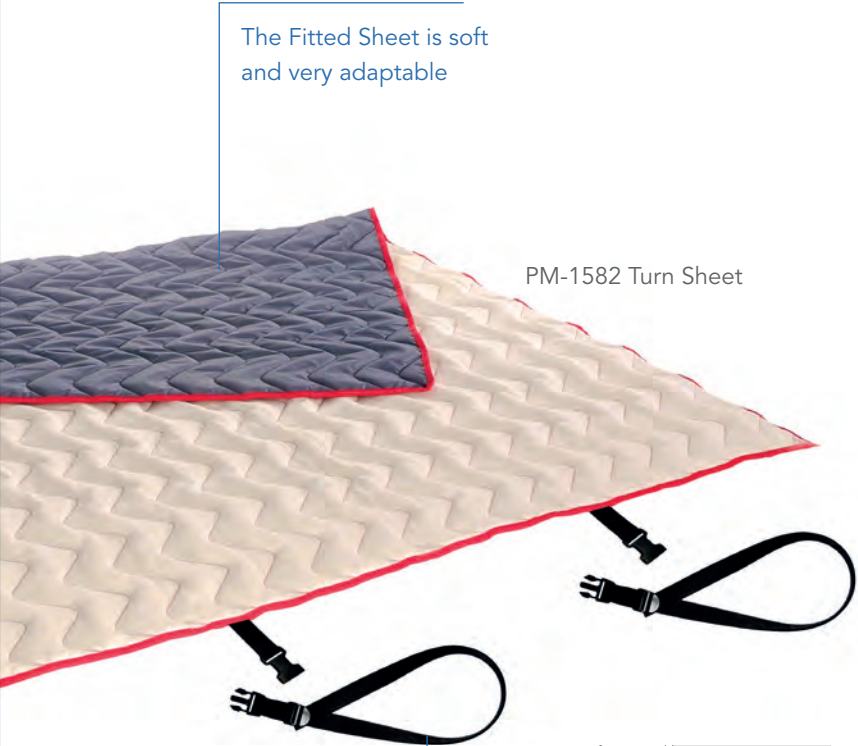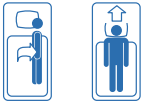

The Fitted Sheet is soft and very adaptable

Three straps enable attachment to the head-board

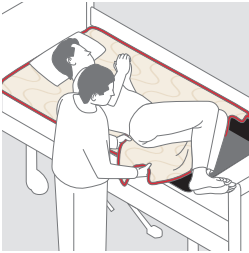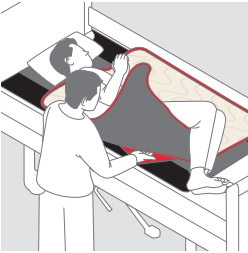

PM-1584 Fitted Sheet



PM-7010 Non-Slip-Mat



PM-8010 Bed Wedges

Turn Sheet Set 135	
Item no.	PM-1580
consisting of:	
Turn Sheet	PM-1582
Fitted Sheet	PM-1584
Non-Slip-Mat	PM-7010
Bed Wedges	PM-8010

Turn Sheet	
Item no.	PM-1582
Dimensions	150 x 200 cm
Material	Cotton/natural ticking

Fitted Sheet	
Item no.	PM-1584
Dimensions	200 x 90 cm

Non-Slip-Mat	
Item no.	PM-7010
Dimensions	35 x 25 cm

Bed Wedges	
Item no.	PM-8010
2 x Bed Wedges	25 x 22 x 13 cm
1 x Bed Wedges	50 x 22 x 13 cm

| Special sizes on request | |


Alpha® Turn Sheet

Set Multi Turn


This Alpha Turn Sheet was specially developed for bedridden patients in need of care. Changing position at short intervals to relieve capillary vessel pressure and constantly repositioning and relieving pressure on areas of the body at risk of pressure sores is essential here. The concealed handles and pull loops on the turning sheet enable carers to position the patient in the desired care position without straining their back. Side turning for basic care is possible without causing pain to the patient, as is positioning. The use of breathable textile products means that the patient can remain on the turning sheet at all times. The stop strips attached to the Fitted Sheet and the Non-Slip-Mat ensure the necessary safety. Thanks to the variable use of the positioning wedges, positioning and relief can be changed with minimal effort, enabling targeted and purposeful positioning.



Pull loops for lifting the patient in bed



PM-2098 Turn Sheet

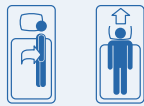

Turn Sheet Set Multi Turn	
Item no.	PM-1695
consisting of:	
Turn Sheet with grips	PM-2098
Fixing Stripes-set	PM-2091
Fitted Sheet	PM-1584
Non-Slip-Mat	PM-7010
Bed Wedges	PM-8010

Turn Sheet with grips		
Item no.	PM-2098	
Dimensions	200 x 100 cm	
Material	Cotton/natural ticking	

Fixing Stripes-set		
Item no.	PM-2091	

Fitted Sheet		
Item no.	PM-1584	
Dimensions	200 x 90 cm	

Non-Slip-Mat		
Item no.	PM-7010	
Dimensions	35 x 25 cm	

Bed Wedges		
Item no.	PM-8010	
2 x Bed Wedges	25 x 22 x 13 cm	
1 x Bed Wedges	50 x 22 x 13 cm	

Special sizes on request
Not suitable for 135° positioning!
(Alternative: Turn Sheets PM-1580)



PM-1584 Fitted Sheet




PM-7010 Non-Slip-Mat




PM-8010 Bed Wedges

Alpha® Turn Sheet

Set Air

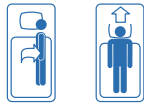

The Alpha Turn Sheet Set Air combines the strengths of the Alpha Turn Sheet Set with the unique possibilities of the Alpha Air. The concealed handles and pull loops on the Turn Sheet enable carers to move patients freely and gently in bed to the desired care position. Thanks to the innovative integration of an air cushion with a variably adjustable air volume, even small movements by the carer or the patient themselves are sufficient to enable them to move to the side.

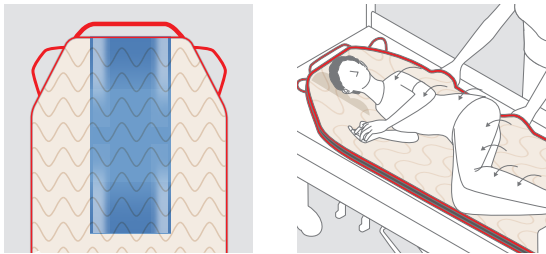

PM-2099.1 Turn Sheet

The air cushion is completely removable and has a washable cover.



PM-1584 Fitted Sheet



PM-7010 Non-Slip-Mat



PM-8010 Bed Wedges

Turn Sheet Set Air	
Item no.	PM-1699
consisting of:	
Turn Sheet with pocket	PM-2099.1
Fixing Stripes-set	PM-2091
Fitted Sheet	PM-1584
Air Cushion	PM-1232.1
Cover for Air Cushion	PM-1233
Air Pump	PM-1234
Non-Slip-Mat	PM-7010
Bed Wedges	PM-8010

Turn Sheet with pocket		
Item no.	PM-2099.1	
Dimensions	200 x 100 cm	
Material	Cotton/natur ticking	

Fixing Stripes-set		
Item no.	PM-2091	

Fitted Sheet		
Item no.	PM-1584	
Dimensions	200 x 90 cm	

Air Cushion		
Item no.	PM-1232.1	

Cover for Air Cushion		
Item no.	PM-1233	

Air Pump		
Item no.	PM-1234	

Non-Slip-Mat		
Item no.	PM-7010	
Dimensions	35 x 25 cm	

Bed Wedges		
Item no.	PM-8010	
2 x Bed Wedges	25 x 22 x 13 cm	
1 x Bed Wedges	50 x 22 x 13 cm	



Alpha® Support Belt

The Support Belt is made of padded polyamide material with non-slip properties inside. It is placed round the patient's and/or caregiver's waist and closed with a buckle. The Support Belt has a number of grips so the carer can control the flow of movement, for example when transferring a patient to a chair, and support insecure patients. The belt is also a reliable aid for walking. Unsteady patients, for example patients with circulation disorders, can be held properly without having to grab their arms.



Support Belt XS

Item no.	PM-6012
for waist	45–70 cm

Support Belt S

Item no.	PM-6013
for waist	60–95 cm

Support Belt M


Item no.	PM-6014
for waist	72–110 cm

Support Belt L

Item no.	PM-6015
for waist	100–150 cm

Support Belt XL

Item no.	PM-6016
for waist	140–180 cm
Material	Polyamide

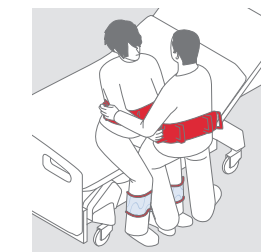

The carer can also wear a belt in order to provide the patient with a good, secure grip

The interior is lined with a non-slip material so it clings to clothing

Vertical and horizontal grips make the handling easier

The buckle is easy to use

The robust outer material is polyamide and prolongs the service life of the belt



Alpha® Raising Aid/Stand-up Support Magic

With the Raising Aid Magic you can securely hold on to and guide the patient when standing up, sitting down and particularly during transfers. When being used, the simple handling can easily be combined with a back-friendly work method and elements from nursing concepts. The Raising Aid Magic can be combined with other aids in various ways, depending on the patients needs and functional requirements. We particularly recommend using the Raising Aid Magic in combination with the leg strap.




Extremely stable carrying handles provide additional safety

With flexible safety harness


Raising Aid Magic with a flexible holding belt	
Item no.	PM-3015
Dimensions	60 x 20 cm
Material	Polyamide
Raising Aid default	
Item no.	PM-3012 (w/o)
Dimensions	60 x 20 cm
Material	Polyamide

With a removable core by means of a clamping fastener and a flexible holding belt

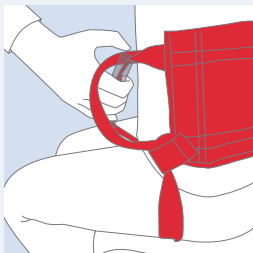

Alpha® Looped Belt

The safety Looped Belt can be used optionally for patients with limited body tension. It enables also especially the safe transfer for patients that sit in very narrow wheelchairs where the flexible holding belt of the Raising Aid Magic cannot be applied to the patient from behind.

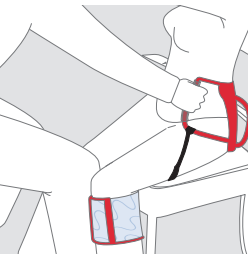

Item no.	PM-3010
Dimensions	110 x 4 cm
Material	Polyester with 2 loops including buckles



With stable buckles



Alpha® Raising Aid-Set Magic & Leg Belt bundle

With the Alpha Raising Aid Magic, you can safely hold and guide people who need assistance when standing up, sitting down and, in particular, when transferring. Its ease of use makes it ideal for combining with back-friendly working methods and elements from care concepts. The Raising Aid Magic with flexible support strap can be combined with other aids as required, depending on the patient's needs and functional abilities. We particularly recommend using the Magic Raising Aid in combination with the leg strap.




Optimally adjusted by extra wide Velcro fastener





Raising Aid-Set	
Item no.	PM-3025
consisting of:	
Raising Aid Magic	PM-3015
Leg Belt	PM-3005

Leg Belt	
Item no.	PM-3005
Dimensions	100 x 22 cm
Material	PES/Carbon

Leg Belt Extension	
Item no.	PM-3007 (w/o)
Dimensions	50 x 22 cm
Material	PES/Carbon

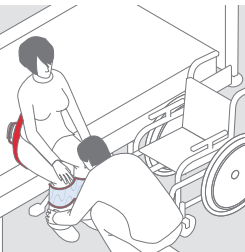


For walking as well the Leg Belt provides safe assistance

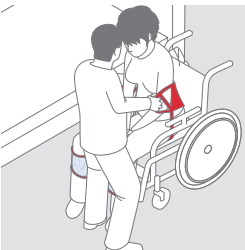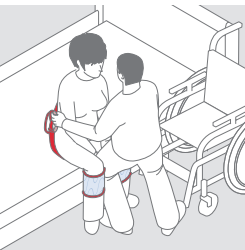

Alpha® Lifting Mat

The Alpha Lifting Mat can be used to move a person in need of care who is lying in bed to the head end and onto their side, enabling transfer from the bed to the shower couch, for example. A durable sliding surface is sewn under the Alpha Lifting Mat. The sliding surface significantly reduces friction between the patient and the transfer surface. The side wings can be easily folded in together with the handles, allowing the Lifting Mat to be positioned optimally under the patient by simply placing it over the semi-lateral position. The inner part of the Lifting Mat can be removed for washing or disinfection. Lifting patients lying on the floor with the Alpha Lifting Mat by several people is safe for the patient and easy on the back for the carers.

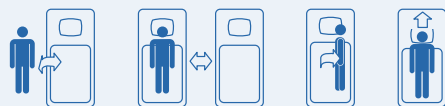

Lifting Mat with 4 grips


Item no.	PM-4012
Dimensions	110 x 56 cm
Material	Polyamide

Lifting Mat with 6 grips

Item no.	PM-4016 (w/o)
Dimensions	150 x 56 cm
Material	Polyamide



The Lifting Mat is available with 4 or 6 grips

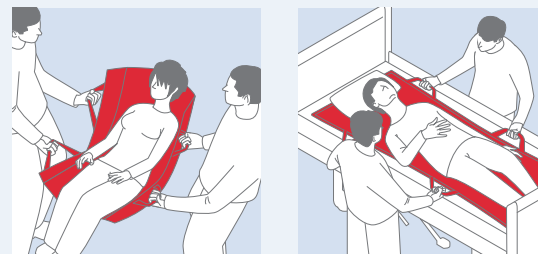

The Lifting Mat not only makes it possible to transfer the patient from a patient support to a bed, it also assists carers when lifting the patient off the floor


Soft padding offers the patient a high level of comfort



The Lifting Mat is foldable for an easier placing into the bed


Alpha® Transfer Seat Swing

The Swing Transfer Seat is ideal for controlled transfers of persons on stairs and other transfer situations, particularly where there is no room for mechanical aids. The weight of the patient can be distributed over up to five persons, even of different heights, to save their backs.

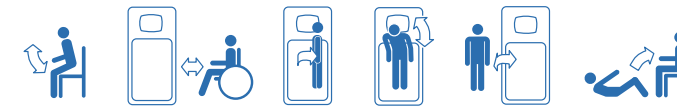


The stable safety handles make the patient transfer safe and easy

Item no.	PM-4035
Dimensions	T 38 x B 43 x H 55 cm
Material	Nylon

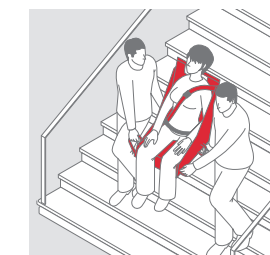



The transfer seat can be used in difficult situations, for example with stairs

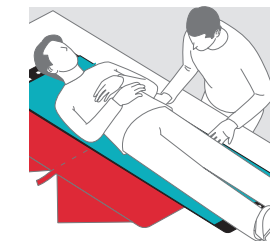


In combination with the Rollboard Vision, a supine patient can be transferred to the opened transfer seat

Alpha® Raising Aid Quick

The Raising Aid Quick is a universal aid with many potential applications. It is ideal for use in conjunction with the sliding board when transferring the patient, e. g. to a wheelchair.

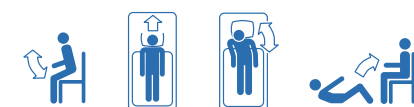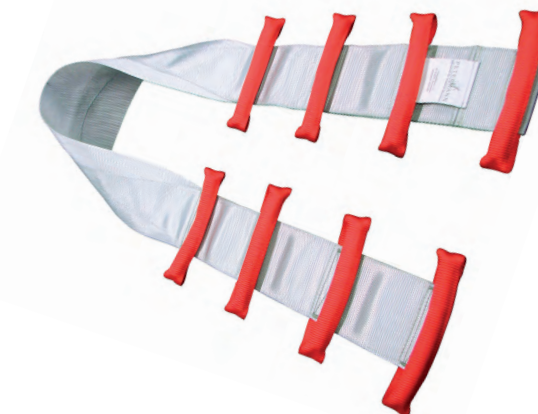


Item no.	PM-6080
Dimensions	145 x 7-13 cm
Material	Polyester

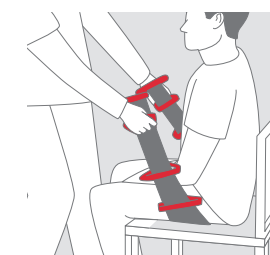



Alpha® Transfer Aid Merlin

The Petermann Merlin Transfer Aid, which can optionally be used with a support and/or looped belt, is a simple and safe aid for assisting people with restricted mobility to stand up and change position. The Merlin Transfer Aid enables carers to transfer patients in an easier, energy-saving and back-friendly manner and it can also mobilise the remaining strength in the patient's arms and legs, if applicable, and also give them a feeling of security. The ingenious design in conjunction with the smooth-running rollers allows a particularly safe and easy handling of the Merlin transfer aid. The use of high-quality rust-free stainless steel also allows disinfection processes even after multiple uses without problems.

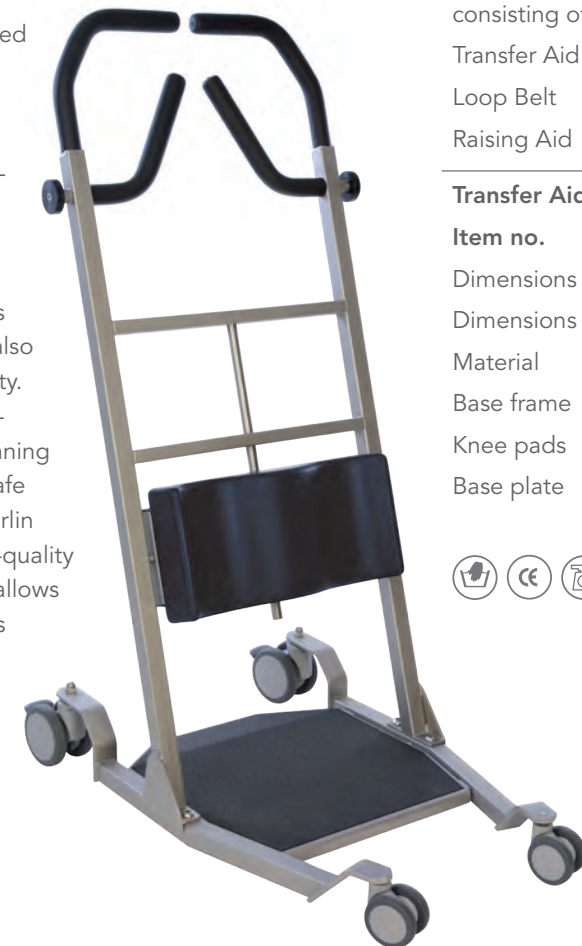

Transfer Aid Merlin, 3-parts	
Item no.	PMM-1000
consisting of:	
Transfer Aid Merlin	PMM-1100
Loop Belt	PM-3010
Raising Aid	PM-6160

Transfer Aid Merlin	
Item no.	PMM-1100
Dimensions Base frame	52,9 x 81 cm
Dimensions Footplate	41 x 46,5 cm
Material	
Base frame	Stainless Steel
Knee pads	PU foam
Base plate	structured plastic film

= 200 kg

Alpha® Support Belt Merlin

The Support Belt Merlin assists with the support of a patient for sitting and standing or for moving a patient into a wheelchair, onto the toilet or onto a chair. It helps the caregiver to control the movement very effectively and obese patients or unstable or nervous patients can be held more securely. Anxious patients are given a sense of security. The Support Belt Merlin is placed round the patient's hips and locked in position with the assistance of the quick-release fastener, the width of which can be adjusted.




Merlin Raising Aid	
Item no.	PM-6160
for waist	70–160 cm
Material	Polyamide

= 200 kg

Alpha® Flexi Board Oval pair

Oval Boards – an innovative aid characterised by versatility and functionality. A wide range of possible uses make everyday life considerably easier for both carers and people in need of support. Low friction, optimised shape and material thickness, and the resulting balanced elasticity enable effortless gliding and insertion. These properties are particularly important for people with limited mobility.

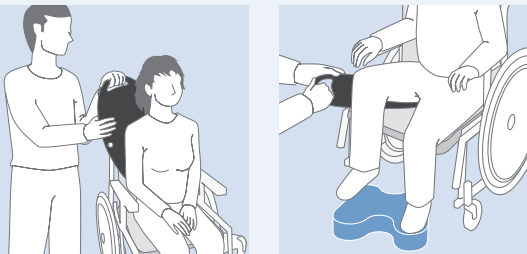

Flexi Board Oval S (pair)	
Item no.	PM-5045
Dimensions 2 x Oval Board	46 x 20 cm
Material	Polyethylene

Flexi Board Oval L (pair)	
Item no.	PM-5040
Dimensions 2 x Oval Board	85 x 40 cm
Material	Polyethylene




Ergonomically shaped handles ensure a good grip

Perfect hanging option thanks to holes on each long side

With the help of the **Flexi Board Oval S** as a set of two, the sitting position in a wheelchair can be corrected flexibly and easily. Transferring with our well-known and proven Alpha Flexible Sliding Boards becomes even easier when the smaller and thinner Flexi Board Oval S is first placed effortlessly under the buttocks of the person being transferred. In a second step, an Alpha Flexible Slide Board is inserted between the Oval Boards. This enables safe transfer, e.g. from the edge of the bed to the wheelchair. In combination with the use of a PM-3005 leg strap, the Flexi Board Oval S provides valuable support when sitting up in bed.

The **Flexi Board Oval L** 2-pack is a unique aid that helps people with sufficient residual strength in their arms to position themselves independently in bed. When using an incontinence pad, the Flexi Board Oval L is slid directly under the pad without having to turn the person in need of care. This makes positioning over the incontinence pad much easier. In combination with an Alpha Rollboard Vision, the Flexi Board Oval L serves perfectly as a pulling aid.

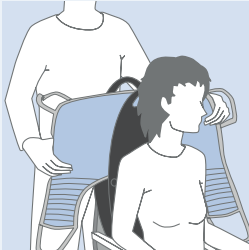

Both sizes of Oval Boards are also ideal for inserting a lifter strap, ensuring safer and more comfortable handling

Alpha® Flexible Sliding Board

The Sliding Board is flexible and has a very smooth surface in order to minimize resistance due to friction. It is used to transfer patients, e. g. from the edge of a bed to a wheelchair, toilet or car seat. The patient does not have to be able to stand up, but can make the change in position while sitting, either with or without assistance.



The anti-slip strips prevent the board from slipping away from the surface

Flexible Sliding Board with rounded corners	
Item no.	PM-5030
Dimensions	74 x 32 cm

Flexible Sliding Board	
Item no.	PM-5050
Dimensions	74 x 32 cm

Flexible Children's Sliding Board	
Item no.	PM-5020 (w/o fig)
Dimensions	58 x 32 cm
Material	Polyethylene



The Flexible Sliding Board makes it possible to transfer a patient to a wheelchair, toilet or car seat

Alpha® Sliding Mat Wheelchair

With the Sliding Mat Wheelchair the user is able to move back in a wheelchair or normal chair on their own or with help from a second person.



Item no.	PM-2070
Dimensions	50 x 45 cm
Material	Cotton/natural ticking



The special material of the Sliding Mat Wheelchair prevents the patient from slipping forward




Alpha® Sliding Mat Wheelchair

The two loops allow a return to the seating position in a particularly comfortable and ergonomic way with the PM-2078 Sliding Mat Wheelchair.



Item no.	PM-2078
Dimensions	135 x 43 cm
Material	Cotton/natural ticking

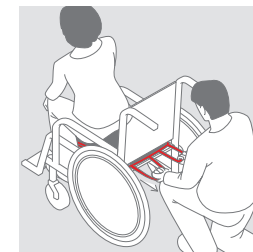

Alpha® Sliding Mat Wheelchair with loops



Item no.	PM-2076
Dimensions	130 x 45 cm
Material	Cotton/natural ticking



With the loops on the Sliding Mat Wheelchair the user can be pulled back by another person without any major effort




The Sliding Mat Wheelchair is attached to the backrest of the wheelchair with loops

The outer surface of the Sliding Mat is made of soft and breathable cotton

Our Website

On our website, you will find informative, helpful application videos and detailed product specifications for numerous products. You can learn more about our corporate values and the high quality standards we set for our products and services. We also provide information about upcoming trade fairs and events organised by Petermann GmbH – an ideal opportunity to get to know us personally and discover our latest innovations.

clear,
informative and
helpful

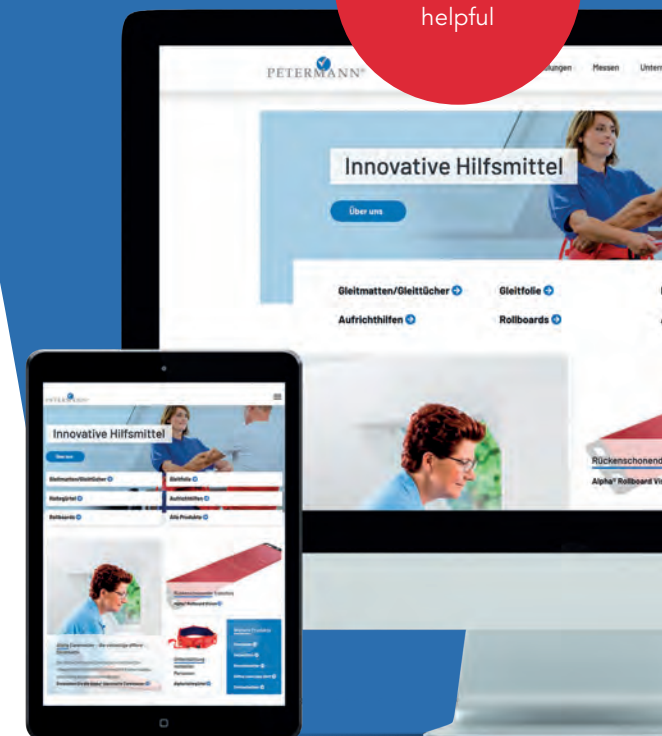

„Visit us at
www.pm-med.com!“




Alpha® Turning Disc

The Alpha Turntable can make it easier for patients who are able to stand to transfer, e.g. between bed and wheelchair. Excellent anti-slip properties enable secure support even on smooth floor coverings. We recommend using the Alpha Turntable in combination with the Alpha Support Belt.

Item no. PM-5080
Diameter 40 cm
Material ABS coated plastic fabric

   = 200 kg

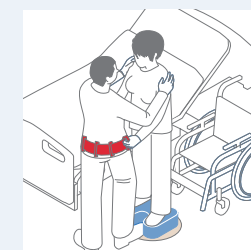

The surface is made of non-slip material in order to prevent the feet from slipping

Load-dependent pivot point for perfect transfer control



The Turning Disc has a grip for carrying it comfortably and spacesaving storage

Alpha® Footstool

The Footstool is used to assist the patient when using the toilet and transferring from bed to wheelchair. It increases the surface of support for the patient. The Footstool can be applied on its own or in combination with the Turning Disc.


Item no. PM-8090
Dimensions 32 x 38 x 8,5 cm
Material Polyethylene

   = max. 200 kg

Protection Cover Footstool

Item no. PM-8091
Material Polyamide

    = max. 200 kg

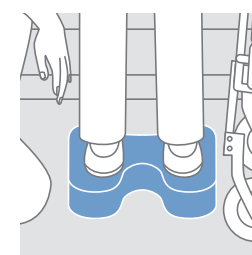

The Footstool helps to overcome height differences, e.g. between bed, wheelchair or toilet


Alpha® Non-Slip Cushion

The Non-Slip Cushion is placed underneath the thighs of a patient seated in bed. A gentle and comfortable supporting and securing function is achieved through the anti-slip properties of the surface. It significantly prevents the patient from sliding into bed from a sitting position.

Item no.	PM-8050
Dimensions	40 x 50 cm
Material	Polyethylene



The Non-Slip Cushion secures the sitting position after the headboard is raised

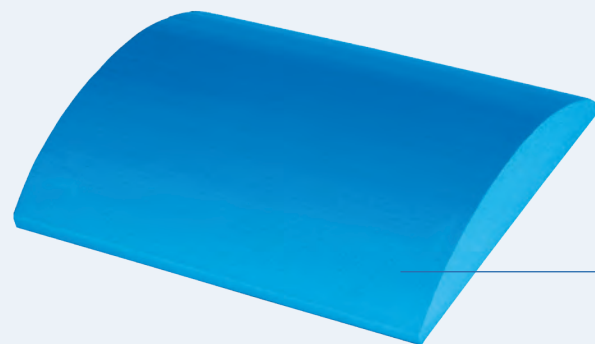


The Non-Slip Cushion prevents a patient seated in bed from slipping down

Non-Slip-Mat

The Non-Slip-Mat prevents the feet from slipping in bed. A Non-Slip-Mat increases the frictional resistance by 100%. When placed underneath the patient's feet in bed, the Non-Slip-Mat helps to prevent slipping while repositioning or raising the patient.

Non-Slip-Mat	
Item no.	PM-7010 red
Item no.	PM-7010-1-blue
Dimensions	35 x 25 cm

Non-Slip-Mat large	
Item no.	PM-7011
Dimensions	45 x 38 cm
Material	Thick film

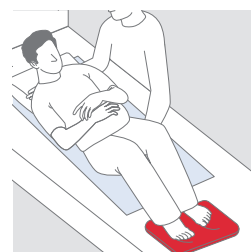

The structure of the mat prevents from slipping


The Non-Slip-Mat makes it possible for a patient moving towards the top end to press with their feet

Alpha® Soft Turning Disc


The Soft Turning Disc is a very simple aid that significantly facilitates lateral movements for seated patients, whether alone or with trusted help. For use in a bed, place the Turning Disc centrally in the bed approx. 25 cm from the edge of the bed. You thus achieve an optimal lying position in the centre of the bed from the turning movement. While in the seating position on the bed always use the height adjustment of the bed to assist with your safety.



The Soft Turning Disc is made of breathable cotton

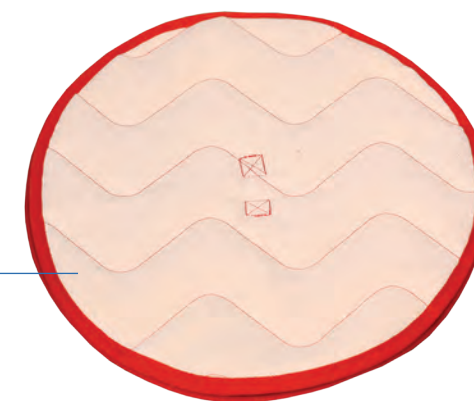

Item no.	PM-5070
Diameter	46 cm
Material	Cotton/natural ticking

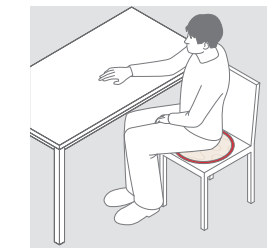


Facilitates sideward rotation from a sitting position

Alpha® Bed Ladder

The Bed Ladder acts as a raising aid and it is made of polyester exclusively. Unlike standard bed ladders, the Petermann Bed Ladder does therefore not lie heavily on the bed and it can be put down quietly. It supports the patient with a physiological pattern of movement when raising the upper part of the body and finally when reaching an upright seating position.

The Bed Ladder enables patients to independently pull themselves into a sitting position in bed

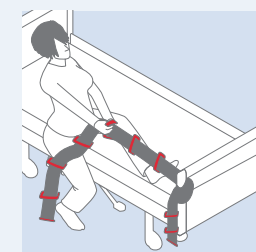

The padded handles allow a secure hold. The grip loops are staggered to make it easier for use

Item no.	PM-6090
Dimensions	170 x 8 cm
Material	Polyester



The Bed Ladder is attached to the bed with a loop

Alpha® X-Ray Mat

The Alpha X-Ray Mat has been specially designed as an X-ray table cover for positioning patients. The comfortable surface and padding provide patients with a comfortable support during X-rays. The surface is light and easy to disinfect. Thanks to its high load-bearing capacity, non-slip underside and side handles on the long sides, nursing staff can easily move patients on the X-ray table.



X-Ray Mat 60 cm

Item no.	PMF-8200
Dimensions	190 x 65 cm

X-Ray Mat 85 cm (w/o)


Item no.	PMF-8220
Dimensions	190 x 85 cm
Material	PVC/PES



The handles on the side of the mat made of radio-paque material are extremely suitable for adjusting the patient's position on the X-ray table

Alpha® X-Ray Crawler

The Alpha X-Ray Crawler is the perfect aid for the professional use of the X-ray cassette on the patient. The innovative concept of the rotating X-Ray Crawler facilitates and accelerates sliding in the X-ray cassette and represents the next important step in the application of the X-ray cassette on the patient. In order to meet highest hygienic requirements, the Alpha X-Ray Crawler can be easily wipe disinfected. The X-Ray Crawler set also offers a disposable protection foil for the X-ray cassette to allow a save and fast workflow with the Alpha X-Ray Crawler.

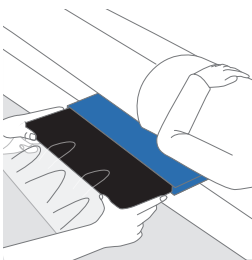

X-Ray Crawler-Set 3-parts

Item no.	PMF-8265
2 x X-Ray Crawler	PMF-8260
1 x Protection Foil	PMF-8261

X-Ray Crawler-Set 6-parts

Item no.	PMF-8266
5 x X-Ray Crawler	PMF-8260
1 x Protection Foil	PMF-8261

X-Ray Crawler

Item no.	PMF-8260
Dimensions	41 x 53.5 cm
Material	Polyamide

X-Ray Crawler-Set XL 3-parts

Item no.	PMF-8365
2 x X-Ray Crawler	PMF-8360
1 x Protection Foil	PMF-8361

X-Ray Crawler-Set XL 6-parts

Item no.	PMF-8366
5 x X-Ray Crawler	PMF-8360
1 x Protection Foil	PMF-8361

X-Ray Crawler XL

Item no.	PMF-8360
Dimensions	47 x 53,5 cm
Material	Polyamide



Protection Foil

Accessories for X-Ray Crawler

Protection foil

Item no.	PMF-8261
Dimensions	45 x 80 cm
Material	LDPE

Protection Foil XL

Item no.	PMF-8361
Dimensions	51 x 80 cm
Material	LDPE

Perforated on roll of 100 pcs

Alpha® Rollboard Vision

The Alpha Rollboard Vision impresses with its confident and back-friendly handling during patient transfers. Uncompromising transfer and care safety is achieved here through a modern combination of hygiene requirements, quality, form and function. All Alpha Rollboards are suitable for use in X-ray and MRI procedures.



Cover machine washable and can also be wipe disinfected



Protection foil



PEM Disposable Cover

Hygienically non-critical and semi-critical applications and places of use
Simple disinfection – cover machine washable up to 60 °C

Hygienically semi-critical and critical applications and places of use
Disposable Protection Foil above the cover for quick, simple handling

See OnlyU® Disposable Protection Foil: Page 53

Hygienically semi-critical and critical applications and places of use
PEM Disposable Cover, can be used without nylon cover, for quick and straightforward handling

See OnlyU® Rollboard Vision PEM: Page 52



Item no.	PMF-9050 Neo Genius
Dimensions	85 x 45 cm, 2 x foldable
Textile Cover	PMF-9051
Item no.	PMF-9075 Ambulance Neo
Dimensions	85 x 45 cm, foldable
Textile Cover	PMF-9076
Item no.	PMF-9085
Dimensions	95 x 50 cm, foldable
Textile Cover	PMF-9086/96
Item no.	PMF-9095
Dimensions	95 x 50 cm, not foldable
Textile Cover	PMF-9086/96
Item no.	PMF-9150
Dimensions	150 x 60 cm, foldable
Textile Cover	PMF-9151
Item no.	PMF-9155
Dimensions	150 x 60 cm, not foldable
Textile Cover	PMF-9151
Item no.	PMF-9170
Dimensions	180 x 50 cm, foldable
Textile Cover	PMF-9171
Item no.	PMF-9180
Dimensions	195 x 50 cm, foldable
Textile Cover	PMF-9181

Item no. PMF-9185
See page 53 for further information

Material Rollboard Polyethylene
Covermaterial Polyamide



Rollboard
Icons: hand, CE, and weight icon = 200/300 kg

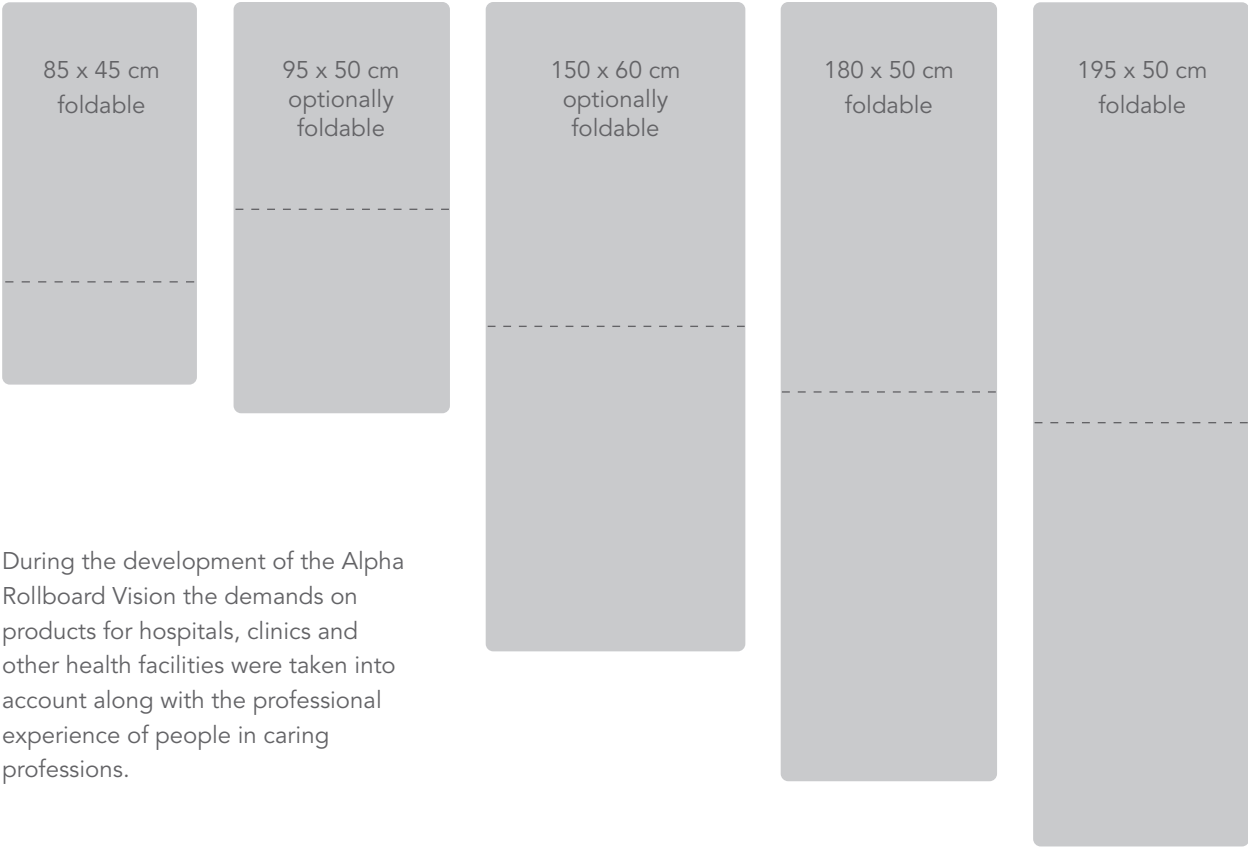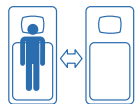


During the development of the Alpha Rollboard Vision the demands on products for hospitals, clinics and other health facilities were taken into account along with the professional experience of people in caring professions.




The patented core of the Rollboard Vision offers numerous advantages; the inner surface enables a quick and thorough disinfection for example

Impacts are absorbed by the flexible Rollboard Vision core; the Rollboard Vision thus provides secure stability and a load-bearing capacity of up to 300 kg, as tested by independent testing bodies

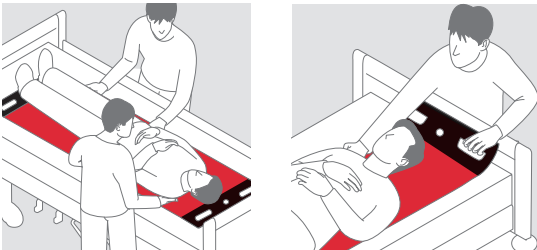

Alpha® Rollboard Vision with cover with handles and grips

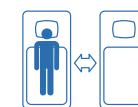

Flexible rollboard cover with
extra large bridging capabilities



Loops enable low contact transfers

Item no. PMF-9185
Dimensions 180 x 50 cm, foldable
with loops and handles


Replacement cover PMF-9186

Material Rollboard Polyethylene
Cover material Nylon

Cover



Rollboard

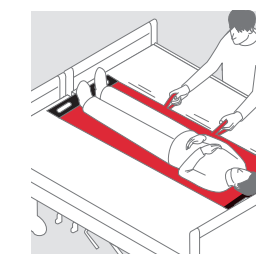


Alpha® Wall Rack for Rollboard Vision

The Wall Rack is designed for easy storage of rollboards. It can be mounted on the wall at a suitable height. Two rollboards can easily be put in the rack and removed for use.



Item no. PMF-0950
Dimensions H 75 x B 45 x T 9 cm
Material Stainless Steel

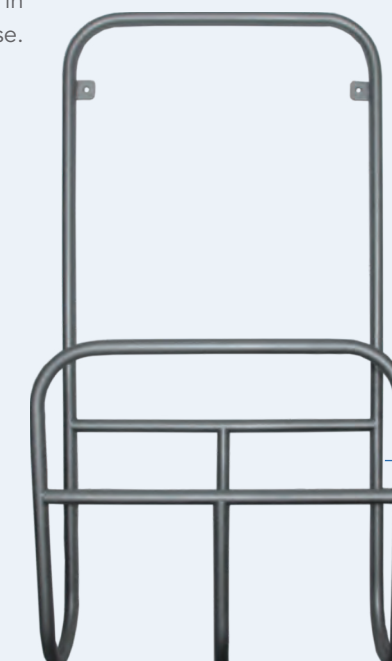

The rollboard can be stored in the
Wall Rack so it is close at hand



OnlyU® Rollboard Vision PEM

Nowhere is the OnlyU Rollboard Vision PEM more impressive than in its confident and back-friendly handling of patient transfers. Uncompromising transfer and care safety is achieved here through a modern combination of hygiene requirements, quality, form and function. All OnlyU Rollboards are suitable for use in X-ray and MRI procedures.



Hygienically non-critical and semi-critical applications and places of use
Simple disinfection – cover machine washable up to 60 °C

See Alpha® Rollboard Vision: Page 46




Hygienically semi-critical and critical applications and places of use
Disposable Protection Foil above the cover for quick, simple handling



Hygienically semi-critical and critical applications and types of use
PEM Disposable Cover, can be used without nylon cover, for quick and straightforward handling

Item no. Set	PMF-9089
Corpus PMF-9085	95 x 50 cm, foldable
PEM Disposable Cover	PM-9090
Item no. Set	PMF-9099
Corpus PMF-9095	95 x 50 cm, not foldable
PEM Disposable Cover	PM-9090
Item no. Set	PMF-9153
Corpus PMF-9150	150 x 60 cm, foldable
PEM Disposable Cover	PMF-9152
Item no. Set	PMF-9157
Corpus PMF-9155	150 x 60 cm, not foldable
PEM Disposable Cover	PMF-9152
Item no. Set	PMF-9175
Corpus PMF-9170	180 x 50 cm, foldable
PEM Disposable Cover	PMF-9172
Item no. Set	PMF-9189
Corpus PMF-9180	195 x 50 cm, foldable
PEM Disposable Cover	PMF-9172
Material Rollboard	Polyethylene
Material PEM Disposable Cover	Polyethylene


Rollboard
= 200/300 kg

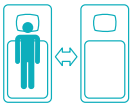

OnlyU® Disposable Protection Foil


Disposable Protection Foil on roll/accessories for Rollboard Vision with textile cover

Item no.	PMF-0191
Dimensions	90 x 54,5 cm
Quantity	110 pcs on roll
Item no.	PMF-0192
Dimensions	180 x 54.5 cm
Quantity	55 pcs on roll

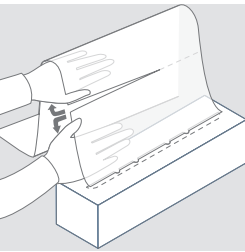


OnlyU® PEM Disposable Cover

Item no.	PMF-9077
Dimensions piece	85 x 47,5 cm
Quantity	50 pcs folded in dispenser bag
suitable for	PMF-9075, PMF-9050
Item no.	PMF-9090
Dimensions piece	90 x 53 cm
Quantity	120 pcs on roll
suitable for	PMF-9085, PMF-9095
Item no.	PMF-9152
Dimensions piece	150 x 61 cm
Quantity	70 pcs on roll
suitable for	PMF-9150, PMF-9155
Item no.	PMF-9172
Dimensions piece	180 x 53 cm
Quantity	60 pcs. on roll
suitable for	PMF-9170
Material	Polyethylene

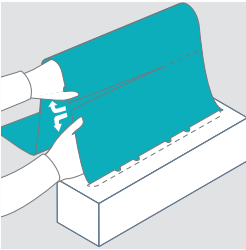


The use of the PEM Disposable Cover reduces the risk of cross-infections and saves on time and costly disinfection processes

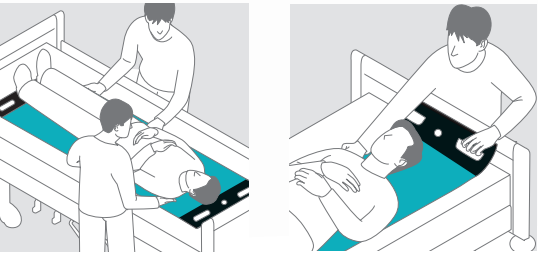



OnlyU® PEM Disposable Sliding Foil

The PEM Disposable Sliding Foil is a tubular sheet which is skin-compatible, tear-resistant, latex free, biocompatible and food-safe. The special composition of the material enhances the high level of slip in daily care. This handling aid is used in clinics and hospitals for all the various patient position changes and transfers. The disposability of the PEM Disposable Sliding Foil means a rapid and particularly hygienic transfer for caregivers so it can be seamlessly integrated into the material administration logistics of a facility. Back-friendly and energy-saving for carers, simple and comfortable for patients – with the PEM Disposable Sliding Foil you reduce friction when you move and/or reposition recumbent or sitting patients.

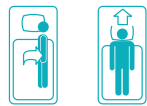

PEM Disposable Sliding Foil 60 cm

Item no.	PM-060.60-100
Dimensions	60 cm x 100 m

**PEM Disposable Sliding Foil 60 cm
incl. dispenser box and foil cutter**

Item no.	PM-060.65-100
Dimensions	60 cm x 100 m


PEM Disposable Sliding Foil 70 cm

Item no.	PM-060.70-100
Dimensions	70 cm x 100 m

**PEM Disposable Sliding Foil 70 cm
incl. dispenser box and foil cutter**


Item no.	PM-060.75-100
Dimensions	70 cm x 100 m

Material	Polyethylene
----------	--------------

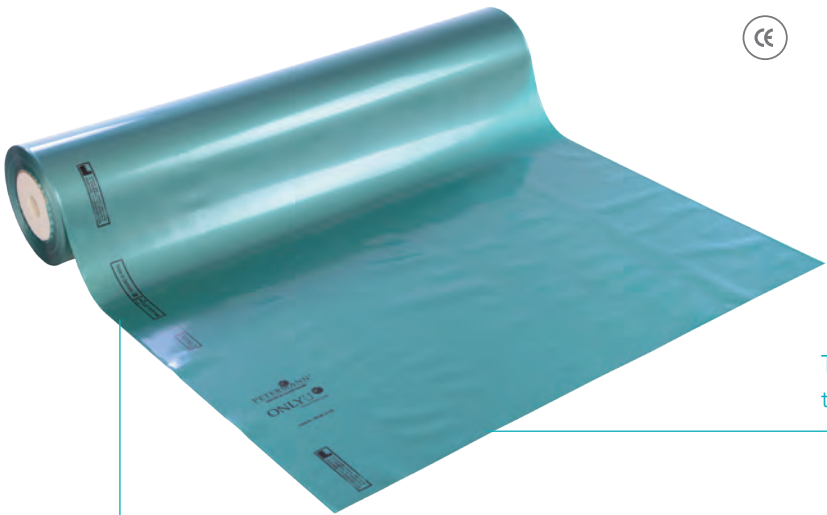

The Sliding Foil can be cut to any desired length

Worked like a tube, the PEM Disposable Sliding Foil can be used just like a Alpha Sliding Mat

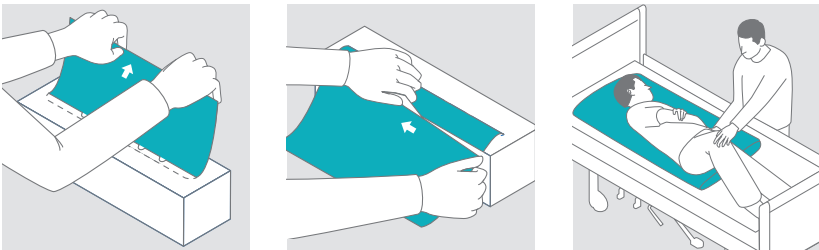


OnlyU® Dispenser for Sliding Foil

Dispenser and cutter for holding a 100 m roll of Sliding Foil



Worked like a tube, the PEM Disposable Sliding Foil can be used just like a Alpha Sliding Mat



Dispenser for roll up to width 70 cm
Item no. PM-060.30

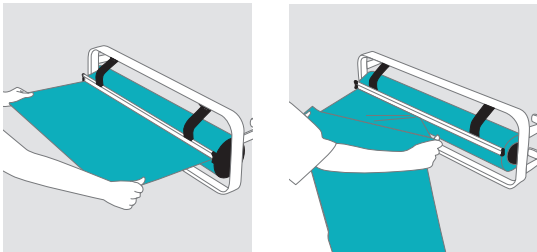

Wall Rack for Foils

The OnlyU Wall Rack was developed for the practical storage of Disposable Protection Foil, PEM Disposable Cover and PEM Disposable Sliding Foil. The Wall Rack can be easily mounted at the appropriate working height. A packaged replacement roll fits into the designated storage compartment.

Item no. PMF-0200.2
Dimensions H 41 x B 80 x T 20 cm
Material Stainless Steel
incl. basket for replacement foil

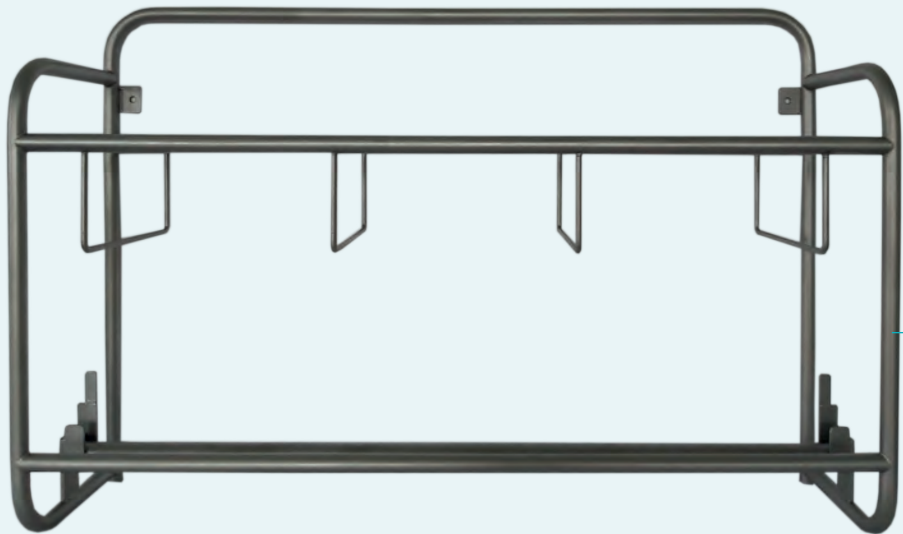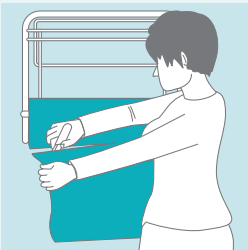

Practical wall rack for foil rolls

Cost Effectiveness

The best offer for you


Modern health care services regard cost effectiveness as the sum of many different factors. In addition to fair purchase prices for aids with long service lives, Petermann GmbH also offers you patient-related disposable products plus many possibilities for calling upon Petermann's know-how in terms of handling products and training.



Blue Motion® Slide

in many sizes

Blue Motion Slides enable patients to be moved in bed. Our Blue Motion Slides also assist with transferring patients to a wheelchair, toilet, chair or shower chair. These can be used in combination with the Alpha Sliding Board, Alpha Magic Raising Aid or Alpha Support Belt.

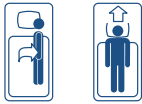


Item no.	BM-0560 *
Dimensions	52 x 60 cm
Item no.	BM-0652 *
Dimensions	60 x 52 cm
Item no.	BM-0752 *
Dimensions	70 x 52 cm
Item no.	BM-S-0790
Dimensions	70 x 90 cm
Item no.	BM-S-0907
Dimensions	90 x 70 cm
Item no.	BM-S-0910
Dimensions	90 x 100 cm
Item no.	BM-S-0911
Dimensions	90 x 110 cm
Item no.	BM-S-1111
Dimensions	110 x 110 cm
Item no.	BM-S-1170
Dimensions	110 x 70 cm
Item no.	BM-S-1190
Dimensions	110 x 90 cm
Item no.	BM-S-1370
Dimensions	130 x 70 cm
Item no.	BM-S-1390
Dimensions	130 x 90 cm
Item no.	BM-S-1970
Dimensions	195 x 70 cm
Material	Polyamide

* recommended for use in wet areas or in combination with the Sliding Board


Special sizes on request



Blue Motion® Slide

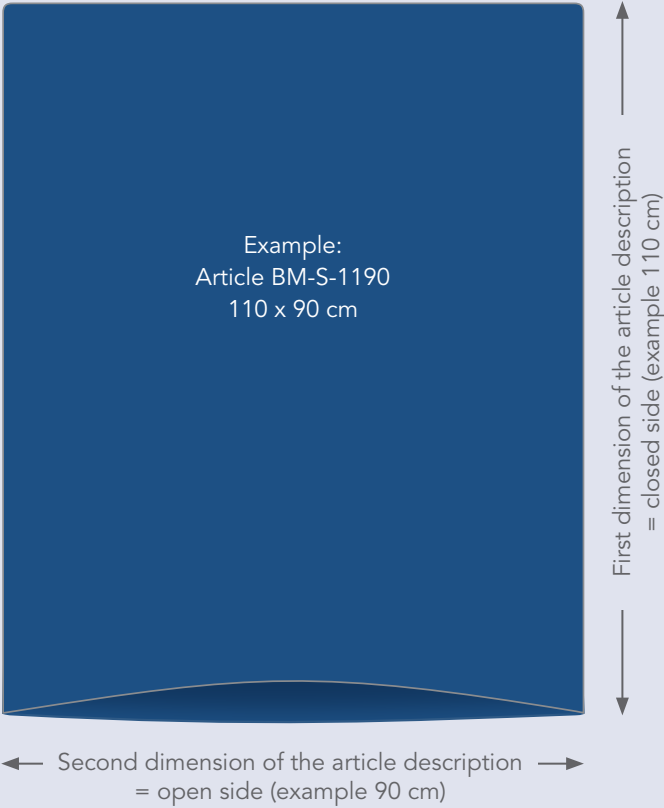

Available
in many sizes

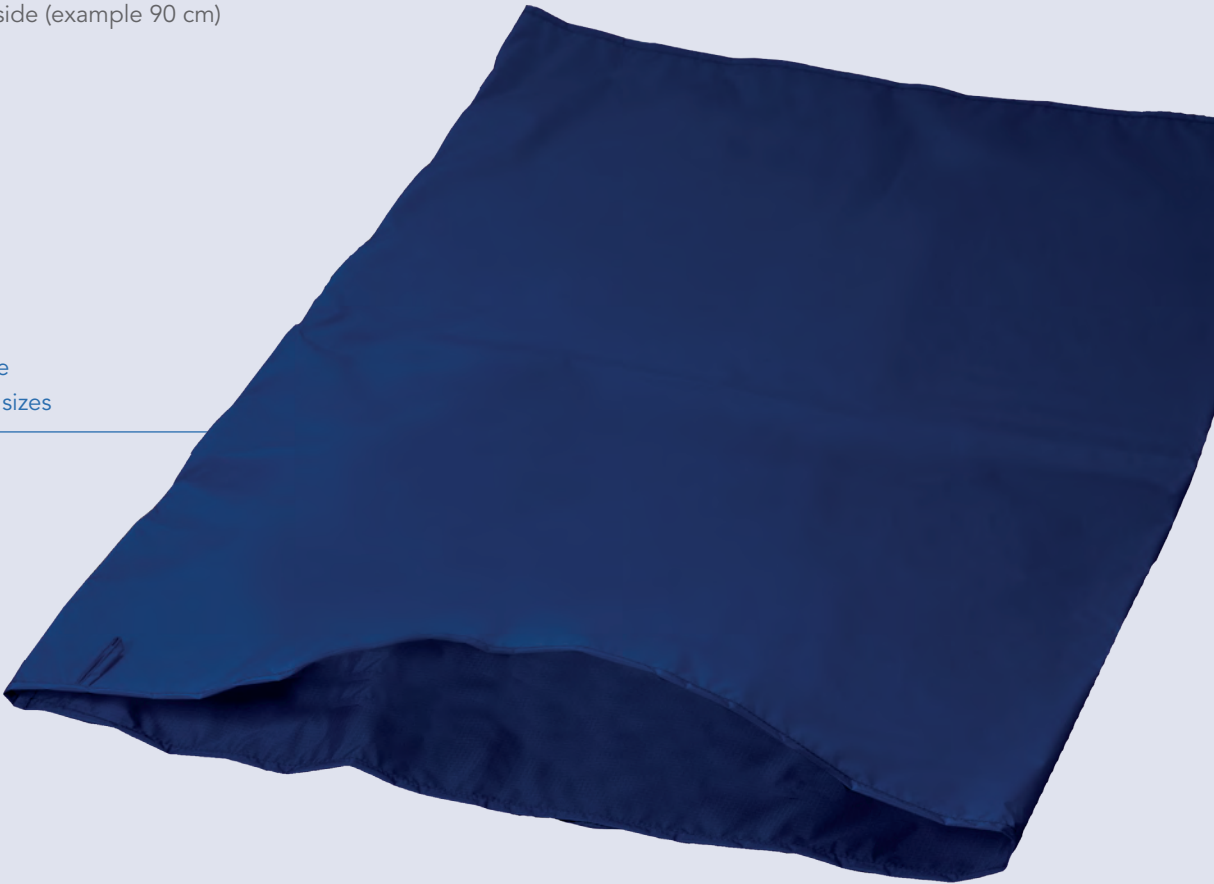

Blue Motion® Non-Slip-Mat

10-pack

The Non-Slip-Mat prevents feet from slipping in bed. The Non-Slip-Mat increases friction resistance by 100 %. Placed under the feet of the person requiring support in bed, the anti-slip mat provides slip resistance during repositioning and lifting.

Item no.	BM-7005
Dimensions	25,5 x 35,5 cm
Material	Thick film



Blue Motion® Sheet

with sewn-in handles

The Blue Motion Sheet does not have a tunnel shape. The sewn-in cord strip allows movement towards the headboard and to the side.

Item no.	BM-S-K-1975
Dimensions	195 x 73 cm
Material	Polyamide
Special sizes on request	



Blue Motion® Sliding Board

The Blue Motion Sliding Board enables the transfer of a person in need of care, e.g. from the edge of the bed to a wheelchair or to the toilet.




Item no.	BM-5000
Dimensions	74 x 24 cm
Material	Polyethylene





EVACUATION

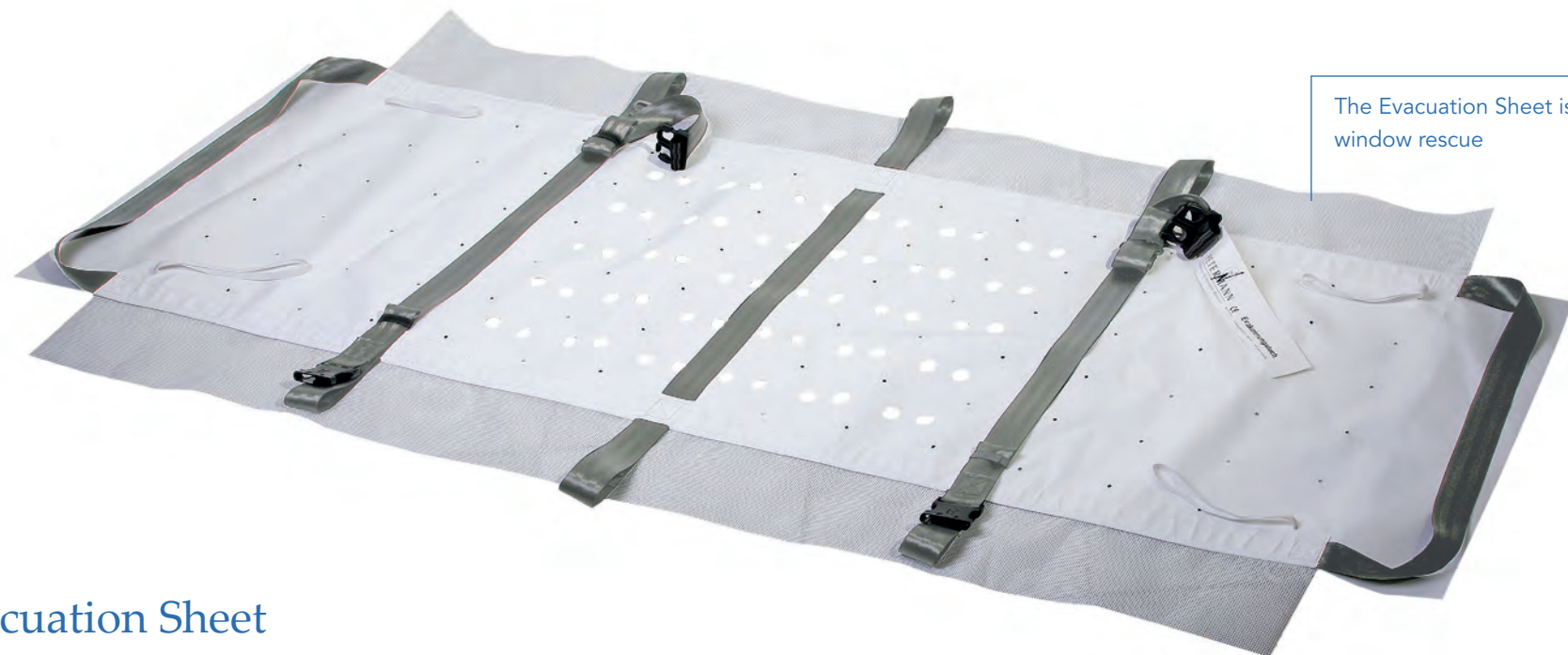

The Evacuation Sheet is ideal for window rescue


Evacuation Sheet

The Evacuation Sheet is used in facilities with patients who are always bedridden. In a rescue situation the Petermann Evacuation Sheet can very easily and safely be dragged over floors, tiles, carpeting, down stairs and out of the hazardous area by a single person. By holding large belt grips at the top or bottom ends and at the centre of the two sides the mattress is pulled out of the bed and across the floor into a safe zone.


Item no.	PM-EVA-7053-H
For bed size	length: 190–205 cm width: 80–100 cm
Item no.	PM-EVA-7053-S
For bed size	length 190–205 cm width: 100–140 cm
Material	PVC coated polyester fabric with ventilation opening

Suitable for all beds/slatted frames/
mattress types

Flame retardant, German TUV approved according to DIN 1865 (Patient handling equipment for rescue services/ambulances), required for Evacuation Sheets

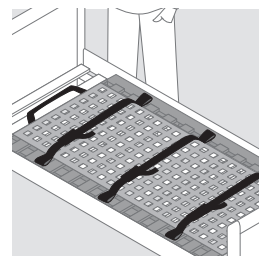


Place the Evacuation Sheet in the middle of the bed directly on the slatted frame.

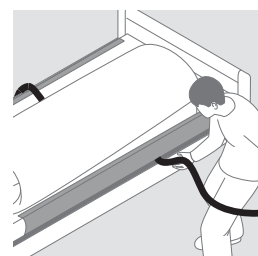


Rescue situation: Pull both plugs of the belts out from under the mattress.



Pull the opposite belt straps out from under the mattress. Connect the opposite belt straps by clicking the buckles into place.

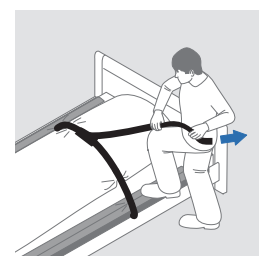


Pull the loose end of each belt as strongly as possible.

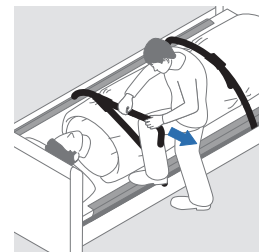


This way the U-shaped mattress protects the person lying inside.



Grab the large handles of the evacuation sheet at head or leg side to pull the sheet together with the mattress from the bed.

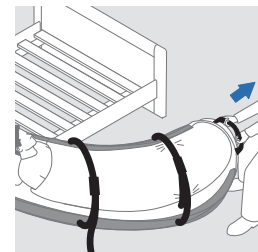

The Evacuation Sheet and mattress slide over the side of the bed and onto the floor.

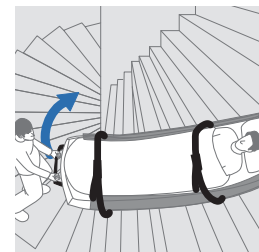

If necessary, you can drag the patient down stairs and into a safe area.



MADE IN GERMANY

Petermann GmbH meets with the requirements
of DIN EN ISO 13485 – medical devices –
Quality management systems. Scope of certificate:
Design, development and distribution of aids
for people of limited mobility

